

UA-767

MAINTENANCE MANUAL

DIGITAL BLOOD PRESSURE MONITOR

AND
A&D Company, Limited




CONTENTS

1. SPECIFICATIONS	2
2. OUTLINE DRAWING	3
3. BLOCK DIAGRAM	4
4. TROUBLESHOOTING	5
TROUBLESHOOTING TABLE	6
5. REPAIR PROCEDURE	7
6. PRESSURE ADJUSTMENT PROCEDURE	8
7. EXPLODED VIEW	10
8. PARTS LIST	11


1. SPECIFICATIONS

1) Measurement Method	Oscillometric
2) Measurement Range	20~280 mmHg (Blood Pressure) 40~200 P/MIN. (Pulse)
3) Accuracy	± 3 mmHg or 2% of measured value (Blood Pressure) $\pm 5\%$ (Pulse)
4) Cuff Inflation	Micropump
5) Cuff deflation	Automatic constant-air-release valve system
6) Rapid Exhaust	Automatic by internal air-release valve
7) Pulse Wave Detection	Manschettor
8) Power source	6VDC, 4 x 1.5V "AA" OR "AM3" batteries
9) Battery life	Approx. 4 months usage of 1 min. per day
10) Weight	Approx. 360 grams.
11) Dimensions	163.7 (W) x 111 (D) x 66.7(H) mm
12) Operating environment	50~100° F. at less than 85% R.H.
13) Storage environment	15~130° F. at less than 95% R.H.
14) Display	Liquid crystal type.

2. OUTLINE DRAWING

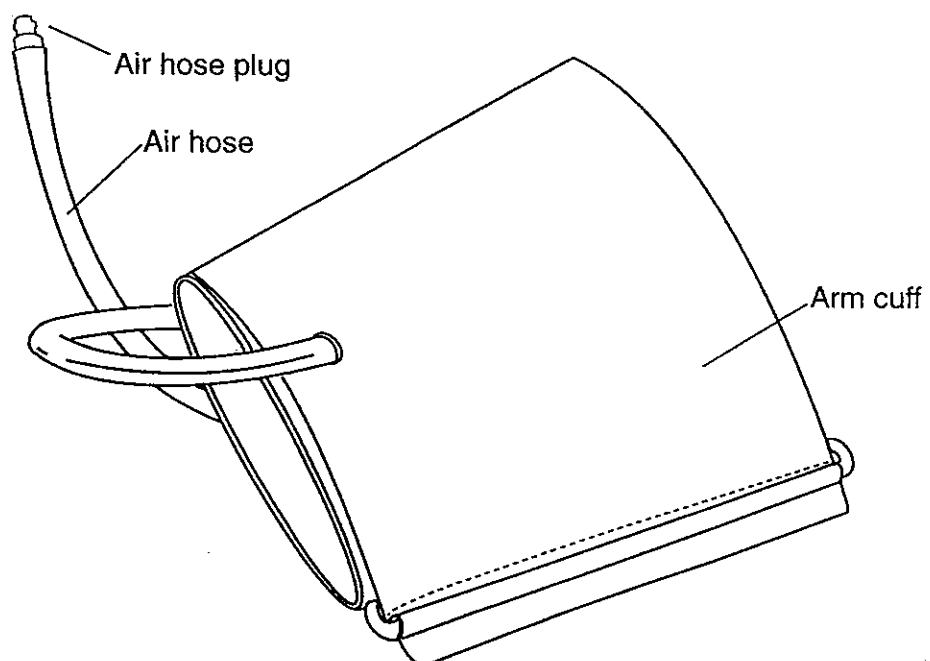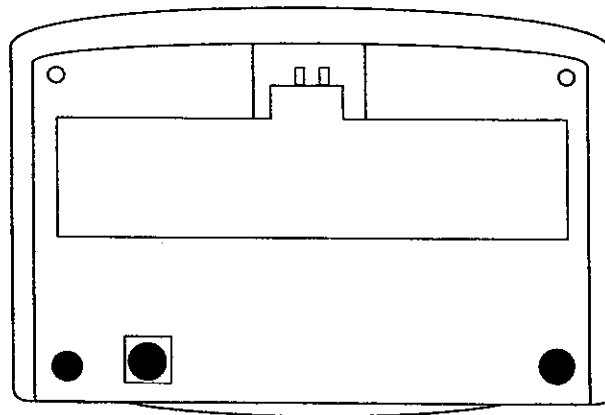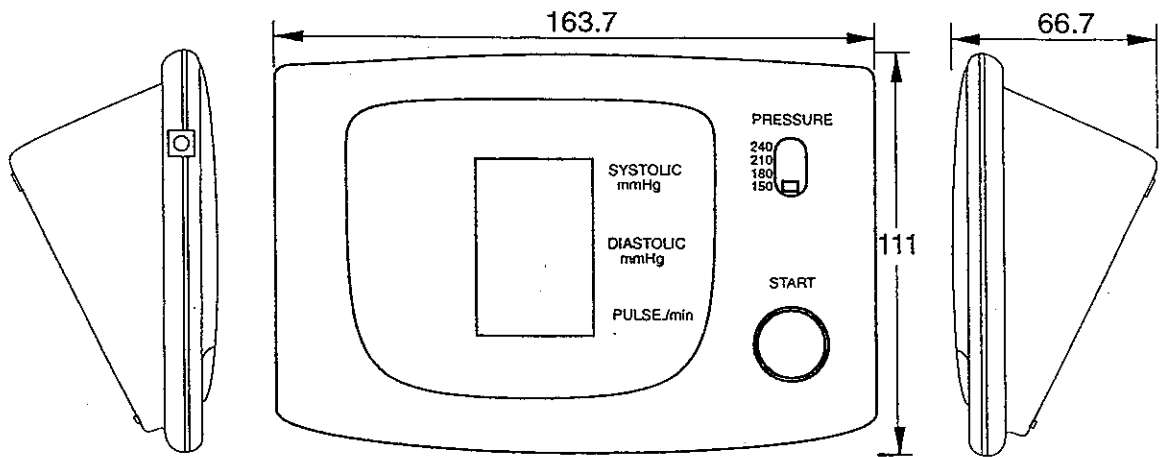

3. BLOCK DIAGRAM

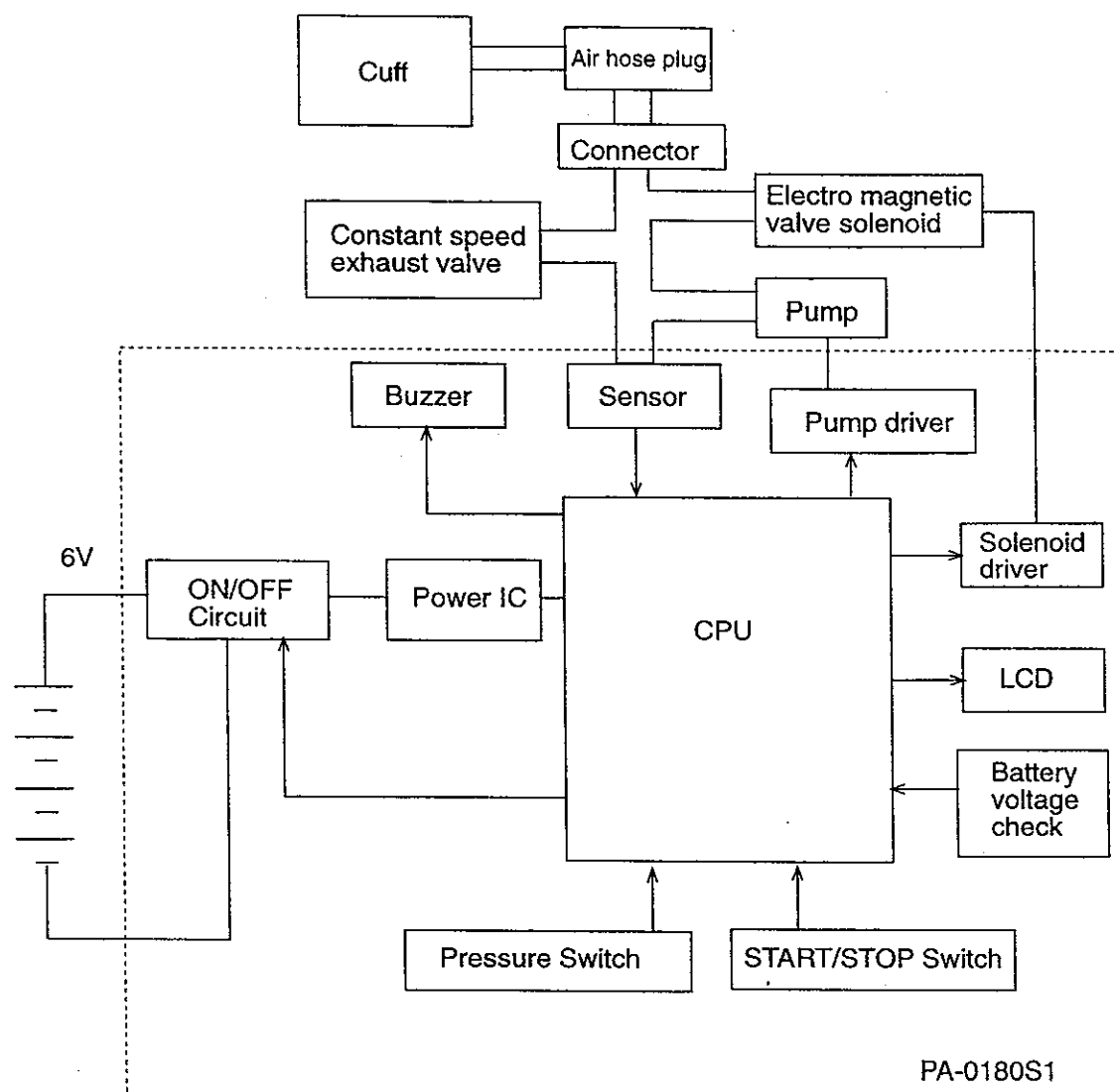


4. TROUBLESHOOTING

This section describes the symptoms, probable causes and solutions to problems. In the case of "can not measure" or "too much error", confirm that the measurement method is correct.

Pressure accuracy should be checked after repair. See "Check Sequence"

Performance check chart

Check the symptoms against the flow chart, find the corresponding number circled on the right side of the chart, then proceed to the troubleshooting table

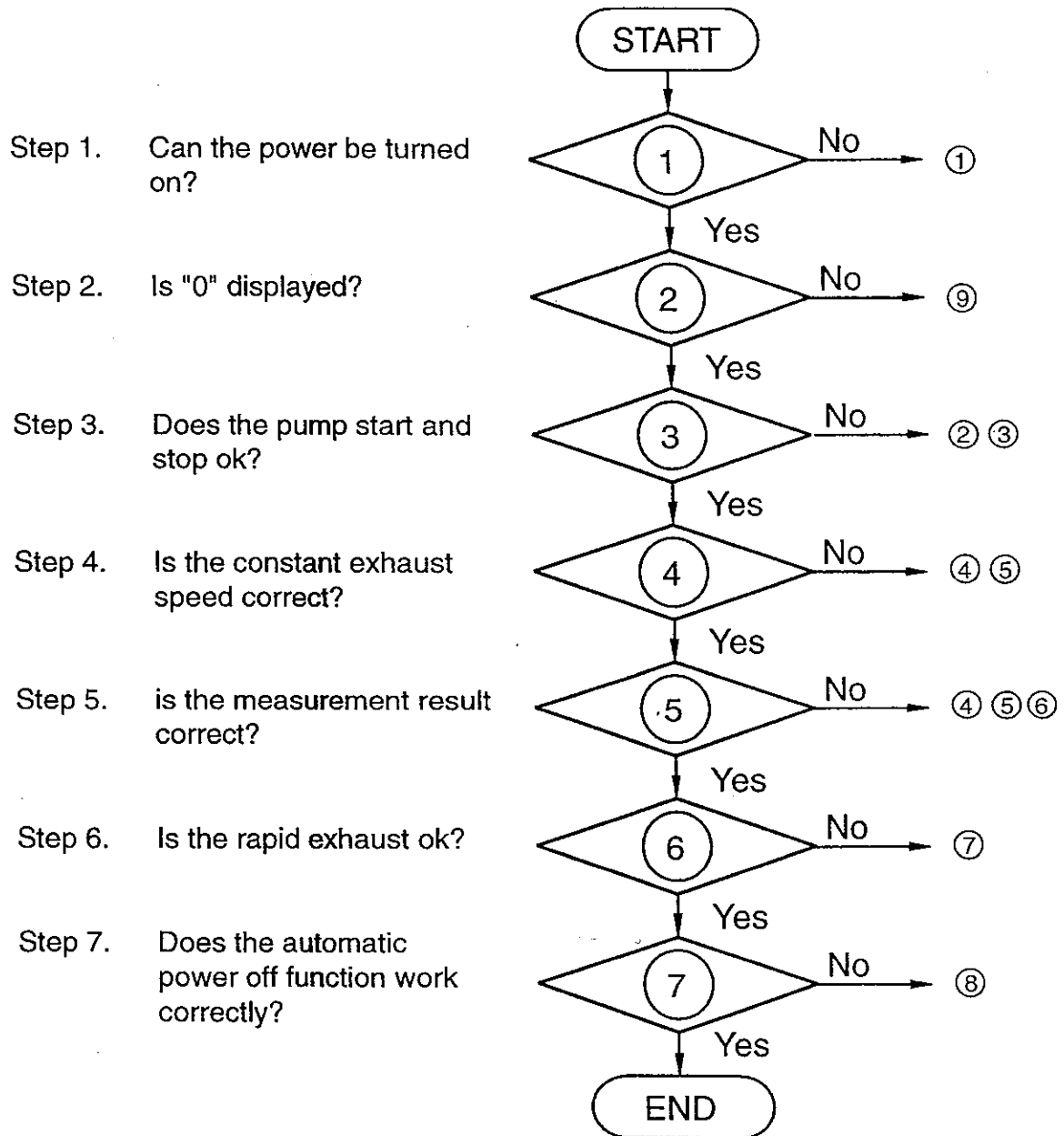


TROUBLESHOOTING TABLE

	Symptom	Probable Cause	Treatment
①	Power does not turn on	Low Battery	Replace battery
		Power lead broken	Resolder leads
		Main board may be defective	Replace main board and adjust pressure reading
②	Pump does not start	Air pump broken	Replace pump
		Connector J1 came off	Reconnect J1 on main board
③	No inflation	Tube came off	Reconnect tubing
		Tubing broken	Replace tubing
		Air connector broken	Replace air connector
		Cuff leaking	Replace cuff
		Constant exhaust valve defective	Replace the constant exhaust valve ass'y.
		Solenoid valve is defective	Replace the Solenoid valve assy.
④	Constant exhaust speed to fast	Constant exhaust valve defective	Replace the constant exhaust valve assy.
		Tubing broken	Replace tubing
		Air connector broken	Replace air connector
⑤	Constant exhaust speed to slow	Constant exhaust valve defective	Replace the constant exhaust valve assy.
		Tubing pinched	Replace tubing
⑥	Pressure reading is incorrect	Pressure reading adjusted incorrectly	Readjust the pressure reading
		Main board may be defective	Replace main board and adjust pressure reading
⑦	Rapid exhaust does not work	Solenoid valve is defective	Replace the Solenoid valve assy.
		Tubing pinched	Replace tubing
⑧	Automatic power off function does not work	Main board defective	Replace main board

5. REPAIR PROCEDURE

Top case removal


- Step 1. Remove the battery compartment cover.
- Step 2. Remove the batteries.
- Step 3. Remove two screws shown in the drawing.
- Step 4. Remove the upper case using caution not to damage the LCD display

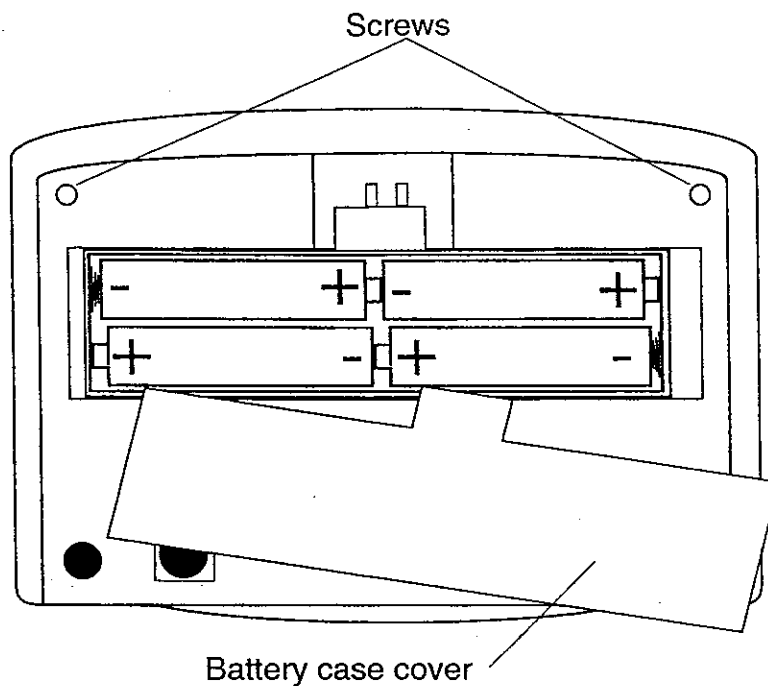

6. PRESSURE ADJUSTMENT PROCEDURE

Test equipment and tools required

- Low capacitance screwdriver
- Manometer or precision mercury pressure gage
- Compressor (rubber ball type)
- T-shaped tube
- Rubber tubing
- Forceps (or hose clamp)

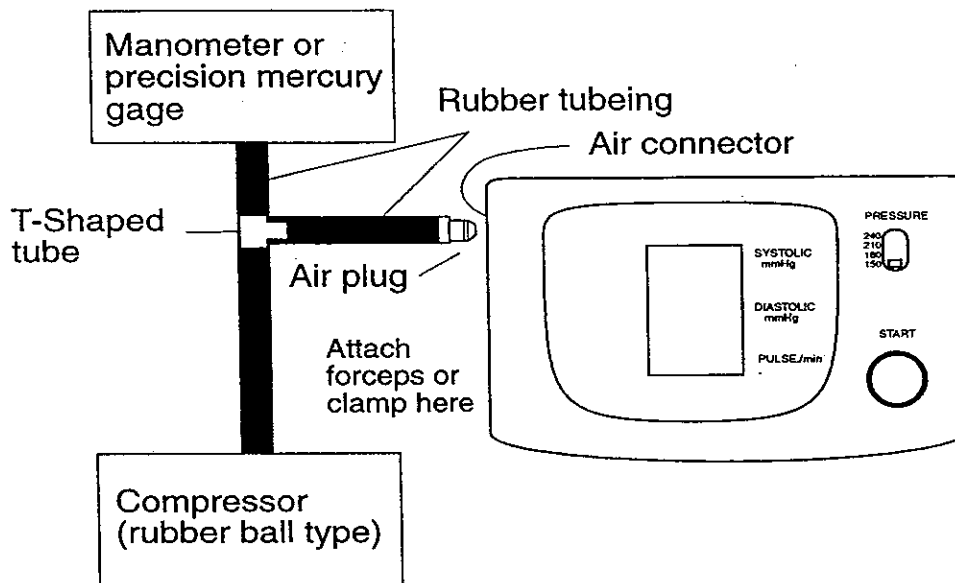

- Step 1. Press the START switch to show all segments of the display. Before "0" is displayed, shift the PRESSURE switch between 150 and 240.
- Step 2. "0" is displayed on the LCD panel
- Step 3. Increase the pressure to 280 mmHg by using the compressor. Use a precision manometer or mercury pressure gage to monitor the air pressure
When using a rubber bulb pump for inflation, close the rubber tube with forceps to maintain the pressure
- Step 4. If the pressure reading is incorrect, remove the top case as described in the repair procedure to provide access to the adjustment screw
- Step 5. Under the following conditions;

Pressure value applied = A mmHg
Reading displayed = B mmHg
Remainder if A - B = C mmHg

Use the formula $A \text{ mmHg} - B \text{ mmHg} = \pm C \text{ mmHg}$

When the result is $+C \text{ mmHg}$, turn the adjustment screw for a reading of $B - C \text{ mmHg}$

When the result is $-C \text{ mmHg}$, turn the adjustment screw for a reading of $B + C \text{ mmHg}$


Example: Pressure value applied = 280 mmHg

Reading displayed = 273 mmHg

$280 \text{ mmHg} - 273 \text{ mmHg} = +7 \text{ mmHg}$

$273 \text{ mmHg} - 7 \text{ mmHg} = 266 \text{ mmHg}$


Adjust the reading to 266 mmHg

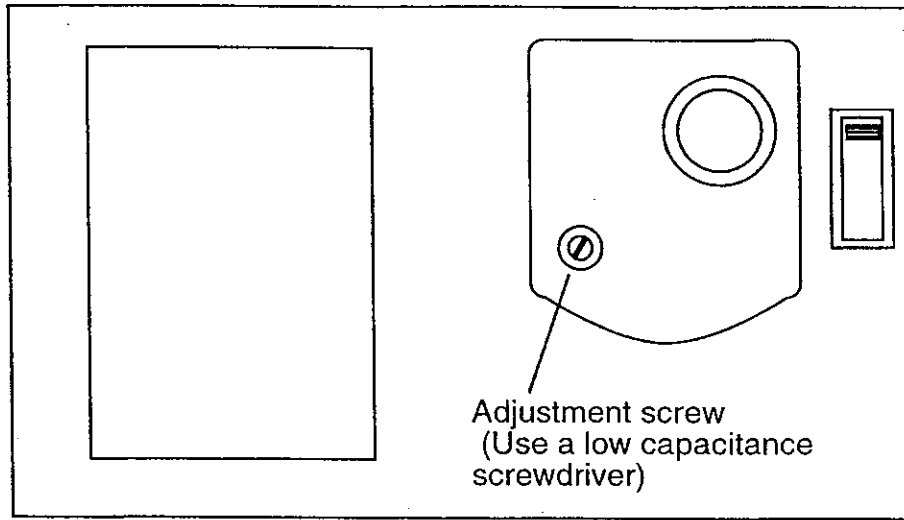

Main board

- Step 6. Reduce the pressure to zero and turn the power off (press the start switch)
- Step 7. Press the start switch turn the power on again (zero is set when the power is first turned on)
- Step 8. Increase the pressure to 280 mmHg by using the compressor and check the accuracy of the display
When the correct reading is obtained, gradually reduce the pressure and confirm that the pressure readings at 150 mmHg and 50 mmHg are within $\pm 3 \text{ mmHg}$
- Step 9. Remove the test setup air plug and attach the cuff to the air connector
Place the cuff on a plastic form about the size of a normal arm
Press the start switch and pressurize the cuff to 160 mmHg
Watch the display reading, the rate of pressure drop between readings should be about 3 mmHg
If necessary adjust the constant exhaust valve for 3 mmHg between readings
- Step 10. Reassemble the case and test the instrument again.

7. EXPLODED VIEW

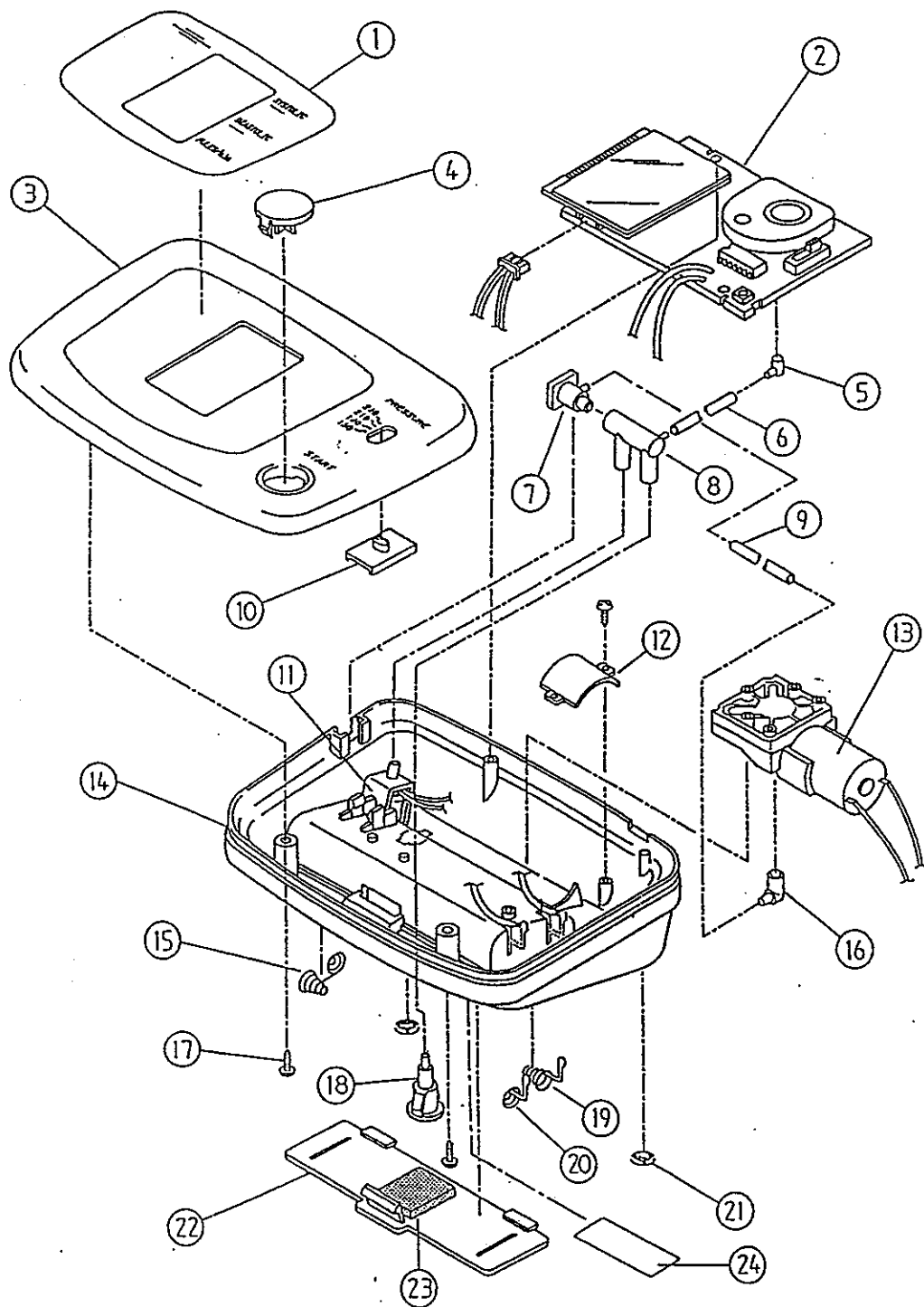

8. PARTS LIST

No.	Parts number	Description	QTY
1	08:4001179	Overlay	1
2	PA:0181S4	PCB Assembly	1
3	07:3000585	Upper case Assembly	1
4	U3-1230	Key cap for Start switch (pink)	1
5	U4-3242-A	Connector for sensor	1
6	TS-23400120TP	Air tube, Ø2.3XØ4.0X120	1
7	U4-5453	Air socket	1
8	UA4-5461	Air connector	1
9	TS-30500070TP	Air tube, Ø3XØ5X70	1
10	U4-5444	Key cap for pressure switch	1
11	LS-TDS-V05B-747	Electro magnetic valve	1
12	U4-5455-A	Pump holder	1
13	UA3-572	Pump Assembly	1
14	U1-194	Lower case Assembly	1
15	U4-1512-A	Battery spring +,-	1
16	U4-3609	L-shape connector for pump	1
17	UZ-0011	Screw M2.3X8mm	3
18	UA4-5348	Constant exhaust valve Assembly	1
19	U4-5472	Battery spring -	1
20	U4-5471	Battery spring +	1
21	U4-5347-A	Rubber foot	2
22	U3-1189	Battery compartment cover	1
23	U4-5342-A	Battery cushion	1
24	08:4002492B	Rating panel	1



A&D Co., LTD.

1-243 Asahi, Kitamoto-shi, Saitama, 364-8585 JAPAN
Telephone: [81] (485) 93-1111 Fax: [81] (485) 93-1119

A&D INSTRUMENTS LTD.

<Authorized Representative Established in the European Community>
Abingdon Science Park, Abingdon, Oxford OX14 3YS ENGLAND
Telephone: [44] (1235) 550420 Fax: [44] (1235) 550485

A&D Medical Pty. Ltd.

32 Dew Street, Thebarton, South Australia 5031
Telephone: [61] (8) 8352 3033 Fax: [61] (8) 8352 7409