

UA-779

MAINTENANCE MANUAL

DIGITAL BLOOD PRESSURE MONITOR

A&D
A&D Company, Limited


CONTENTS

1. SPECIFICATIONS	2
2. OUTLINE DRAWING	3
3. BLOCK DIAGRAM	4
4. TROUBLESHOOTING	5
TROUBLESHOOTING TABLE	6
5. REPAIR PROCEDURE	7
6. PRESSURE ADJUSTMENT PROCEDURE	8
7. EXPLODED VIEW	10
8. PARTS LIST	11


1. SPECIFICATIONS

1) Measurement Method	Oscillometric
2) Measurement Range	20~280 mmHg (Blood Pressure) 40~200 P/min. (Pulse)
3) Accuracy	±3 mmHg or 2% of measured value (Blood Pressure) ±5% (Pulse)
4) Cuff Inflation	Micro pump
5) Cuff deflation	Automatic constant-air-release valve
6) Rapid Exhaust	Automatic by internal air-release valve
7) Pulse Wave Detection	Manschettor
8) Power source	6VDC, 4 x 1.5V "AA" OR "LR06" batteries
9) Battery life	Approx. 4 months usage of 1 measurement. per day
10) Weight	Approx. 320 grams.
11) Dimensions	145 (W) x 130(D) x 56(H) mm
12) Operating environment	+10~40° C. at less than 85% R.H.
13) Storage environment	-20~+70. at less than 95% R.H.
14) Display	Liquid crystal.

2. OUTLINE DRAWING

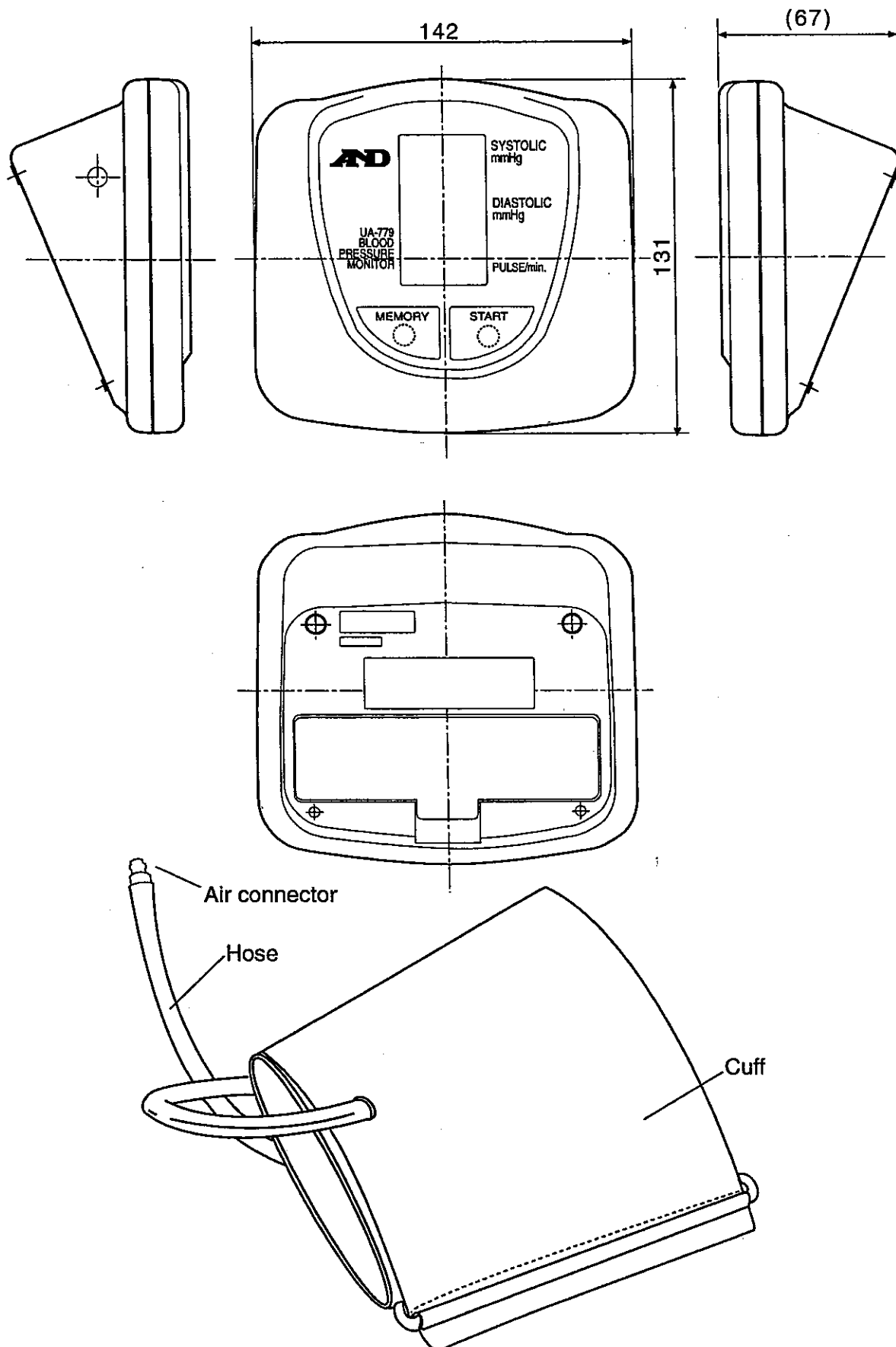

3. BLOCK DIAGRAM

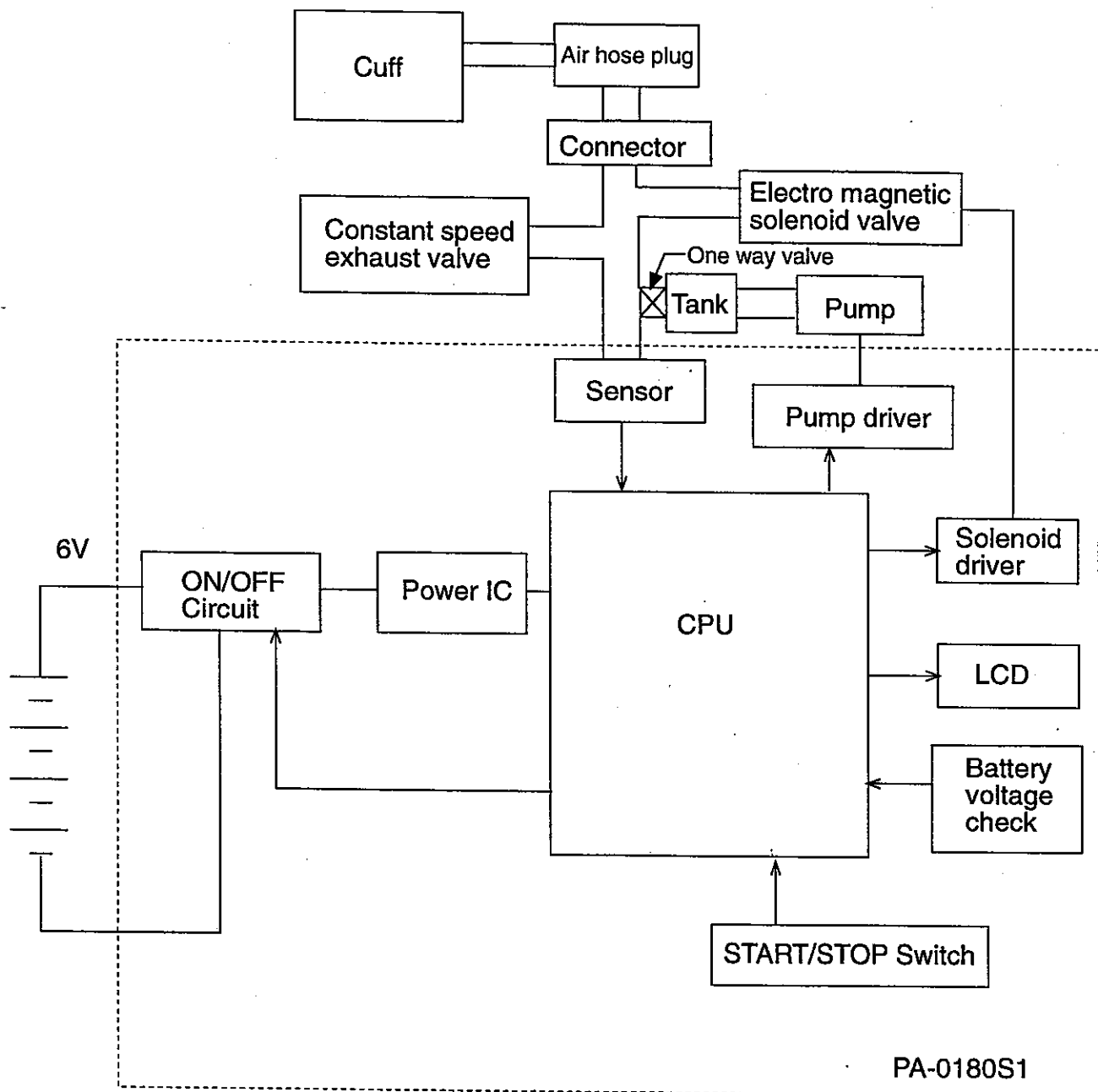

4. TROUBLESHOOTING

This section describes the symptoms, possible causes and solutions to fix THE UNIT.


In the case of "can not measure" or "too much error", check that the measurement method is correct.

Pressure accuracy should be checked after repairing. See "Check Sequence" Performance check chart.


Check the symptoms against the flow chart shown below and find the corresponding number circled on the right side of the chart. Then proceed to the troubleshooting table.

<Check Sequence>


Step 1. Can the power be turned on?

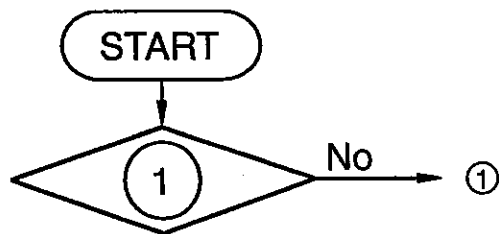


Step 2. Is "0" displayed?

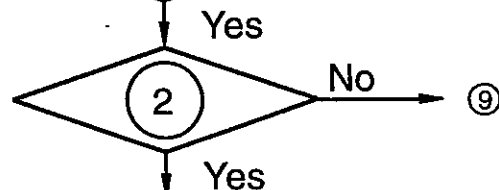


Step 3. Does the pump start and stop ok?

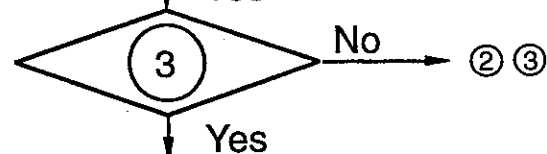


Step 4. Is the constant exhaust speed correct?

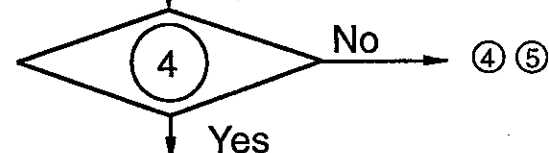


Step 5. Is the measurement result correct?

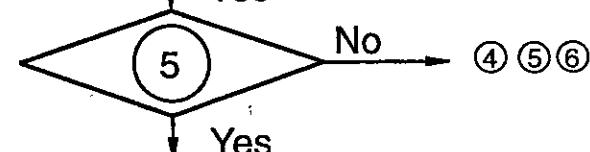

Step 6. Is the rapid exhaust ok?

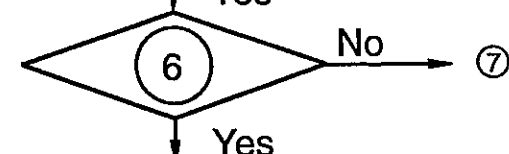

Step 7. Does the automatic power off function work correctly?

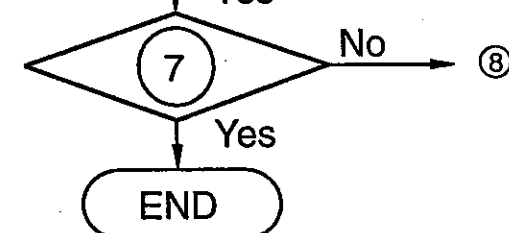


TROUBLESHOOTING TABLE

	Symptom	Probable Cause	Treatment
①	Power does not turn on	Low Battery	Replace the battery
		Power lead broken	Resolder leads
		Main board may be defective	Replace the main board and adjust pressure reading
②	Pump does not start	Air pump broken	Replace the pump
		Connector J1 came off	Reconnect J1 on the main board
③	No inflation	Tube came off	Reconnect tubing
		Tubing broken	Replace tubing
		Air connector broken	Replace air connector
		Cuff leaking	Replace the cuff
		Constant exhaust valve defective	Replace the constant exhaust valve ass'y
		Solenoid valve defective	Replace the Solenoid valve assy.
④	Constant exhaust speed too fast	Constant exhaust valve defective	Replace the constant exhaust valve assy
		Tubing broken	Replace tubing
		Air connector broken	Replace the air connector
⑤	Constant exhaust speed too slow	Constant exhaust valve defective	Replace the constant exhaust valve assy
		Tubing pinched	Replace tubing
⑥	Pressure reading is incorrect	Pressure reading adjusted incorrectly	Readjust the pressure reading
		Main board may be defective	Replace the main board and adjust pressure reading
⑦	Rapid exhaust does not work	Solenoid valve is defective	Replace the Solenoid valve assy.
		Tubing pinched	Replace tubing
⑧	Automatic power off function does not work	Main board defective	Replace the main board
⑨	Pressure sensor unstable	Main board is defective Pressure sensor is defective	Replace the main board

5. REPAIR PROCEDURE

Top case removal

- Step 1. Remove the battery compartment cover.
- Step 2. Remove the batteries.
- Step 3. Remove two screws shown in the battery compartment.
- Step 4. Remove the upper case using caution not to damage the LCD display.

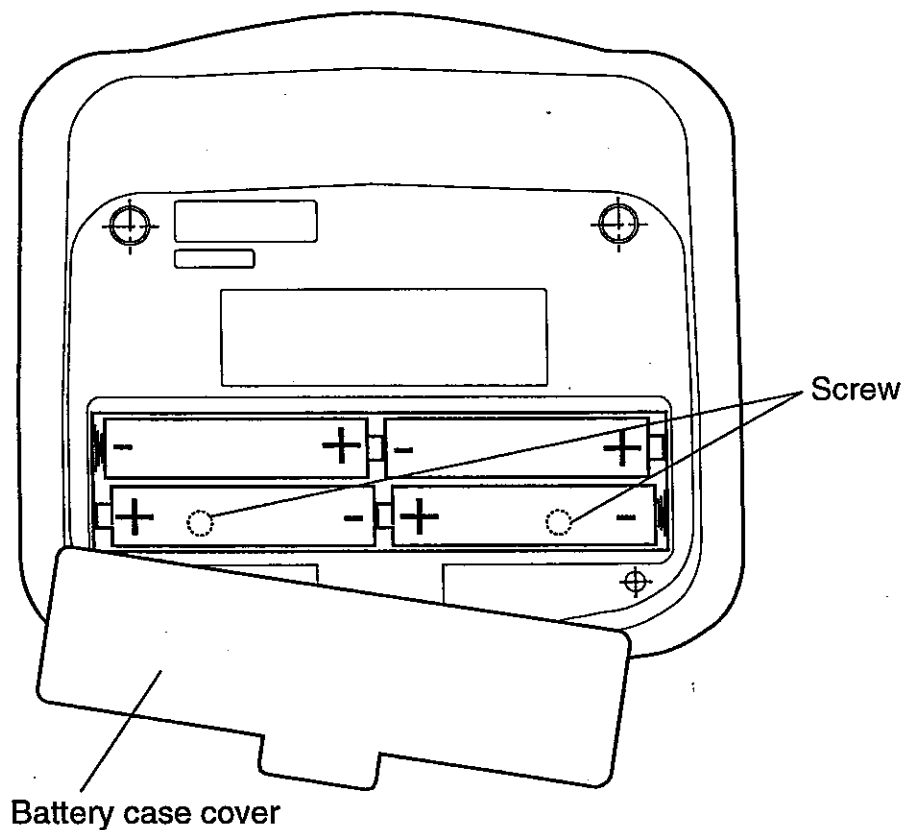


Under side of the monitor

6. PRESSURE ADJUSTMENT PROCEDURE

Test equipment and tools required

- Low capacitance screwdriver
- Manometer or precision mercury pressure gage
- Air compressor (rubber ball type)
- T-shaped tube
- Rubber tubing
- Forceps (or hose clamp)

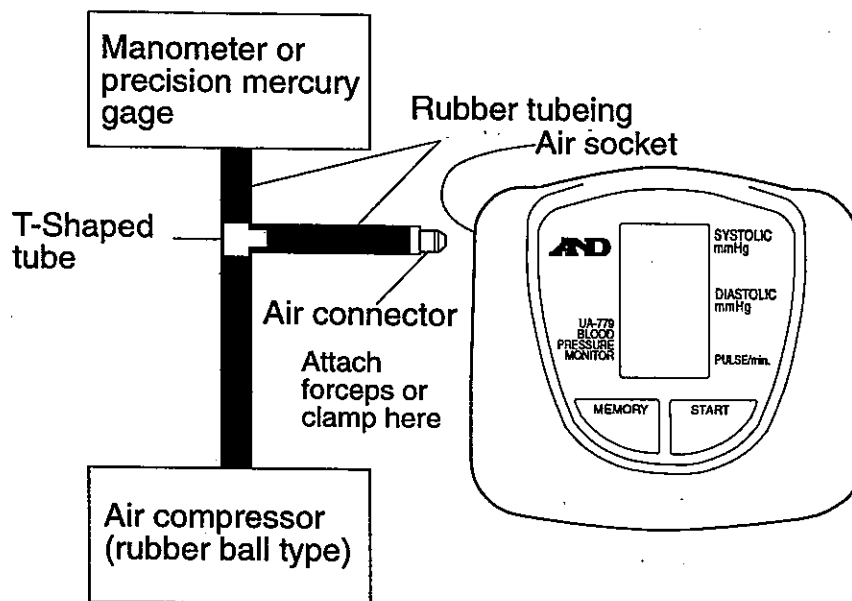

- Step 1. Enter the check mode.
<Entering the check mode>
The battery is inserted while pushing Start key.
- Step 2.

0
0
0

 is displayed on the LCD panel.
- Step 3. Increase the pressure to 280 mmHg by using the compressor. Use a precision manometer or mercury pressure gage to monitor the air pressure.
When using a rubber bulb pump for inflation, close the rubber tube with forceps to maintain the pressure.
- Step 4. If the pressure reading is incorrect, remove the upper case as described in the repair procedure to provide access to the adjustment screw.

Step 5. Under the following conditions;

Pressure value applied = A mmHg

Reading displayed = B mmHg

Remainder if A - B = C mmHg

Use the formula $A \text{ mmHg} - B \text{ mmHg} = \pm C \text{ mmHg}$.

When the reference is +C mmHg, turn the adjustment screw for a reading of B - C mmHg.

When the reference is -C mmHg, turn the adjustment screw for a reading of B + C mmHg.


Example: Pressure value applied = 280 mmHgA

Reading displayed = 273 mmHgB

280 mmHg - 273 mmHg = +7 mmHgC

273 mmHg - 7 mmHg = 266 mmHg

Adjust the reading to 266 mmHg

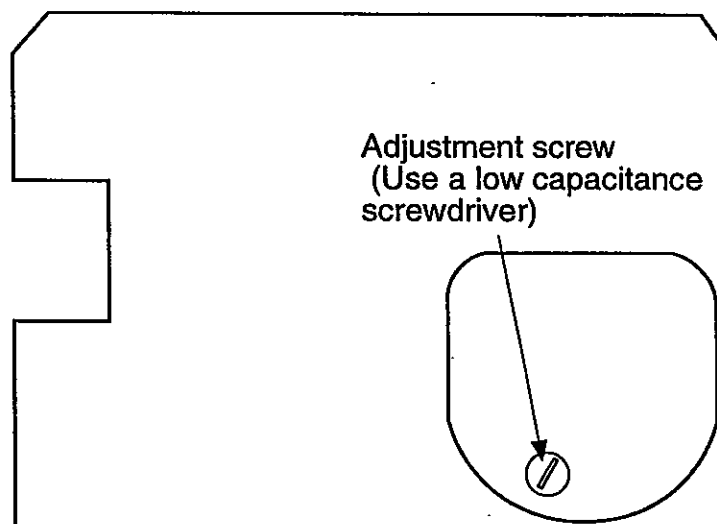

Main board

Step 6. Reduce the pressure to zero and turn the power off. (press the start switch)

Step 7. Enter the check mode. Same procedure Step1.

Step 8. Increase the pressure to 280 mmHg by using the compressor and check the accuracy of the display.

When the correct reading is obtained, gradually reduce the pressure and confirm that the pressure readings at 150 mmHg and 50 mmHg are within ± 3 mmHg against the external gage.


Step 9. Remove the special air plug and attach the cuff to the air connector. Place the cuff on a plastic form about the size of a normal arm. Press the START switch and pressurize the cuff to 160 mmHg by holding.

Watch the display reading. The rate of pressure drop between readings should be about 3 mmHg.

If necessary, adjust the constant exhaust valve for 3 mmHg between readings.

Step 10. Reassemble the case and test the instrument again.

7. EXPLODED VIEW

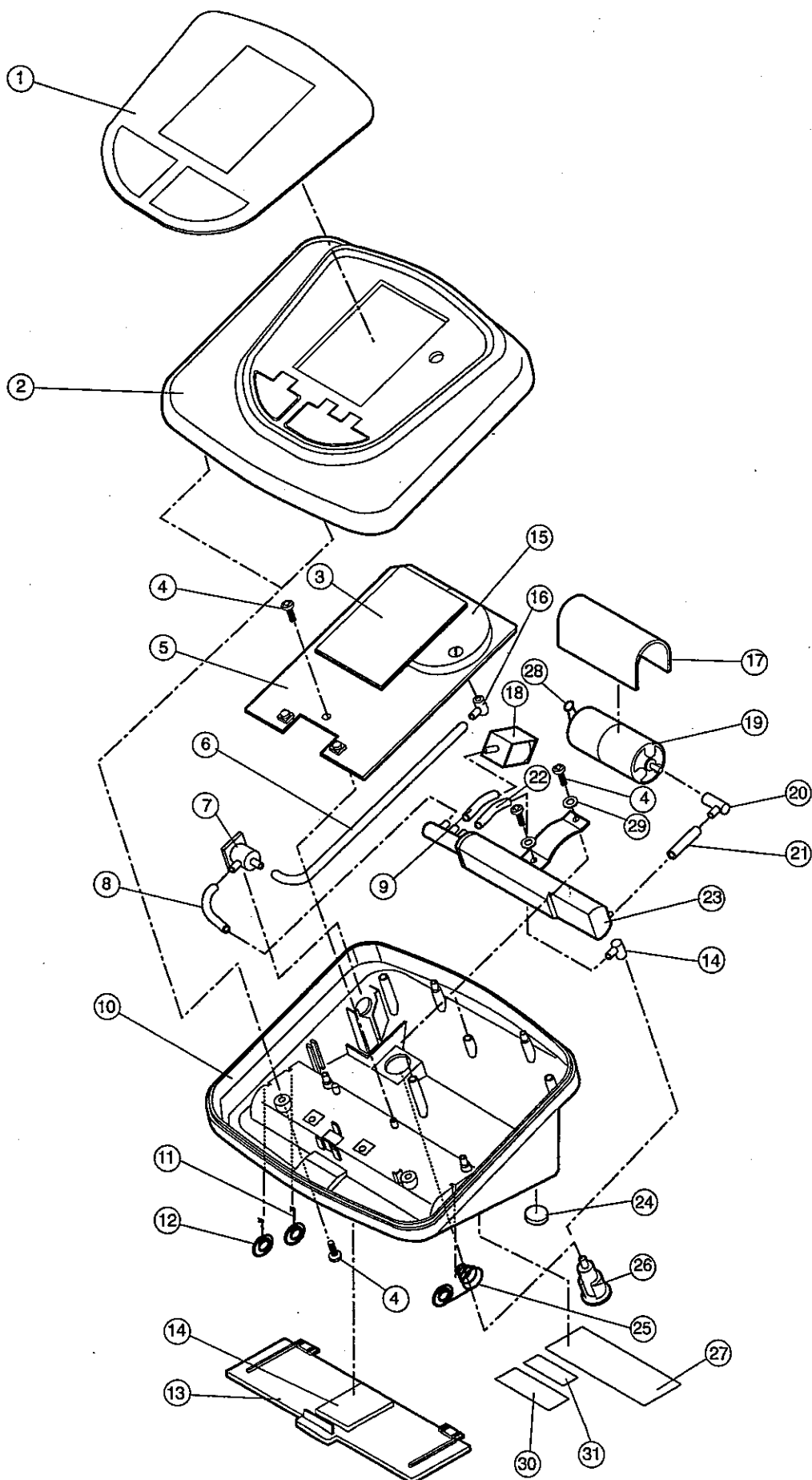

8. PARTS LIST

No.	Parts number	Description	QTY
1	08:4007189	LCD Cover Panel	1
2	07:2000320	Upper Case	1
3	VL-BTJ007-TS	LCD	1
4	UZ4-0011	Screw M2.3 X 8	5
5	PA-0202S1	PCB	1
6	TS-23400130TP	Tube	1
7	U4-4513	Air Socket	1
8	TS-30500025	Tube	1
9	TS-23400017TP	Tube	1
10	07:1000069	Lower Case	1
11	U4-5388	Terminal(+)	1
12	U4-3521	Terminal(-)	2
13	U3-1189	Battery Cover	1
14	U4-5342-A	Battery Cushion	1
15	PA:4000406	Sensor	1
16	U4-3242-A	Sensor Joint	1
17	U4-5346-A	Pump Cushion	1
18	LS-TDS-V05B-764	Relase Valve	1
19	LM-P05H-0007	Pump	1
20	U4-3609	Pump Joint	1
21	TS-30500020TP	Tube	1
22	TS-23400020TP	Tube	1
23	UA4-5353	Tank Assy	1
24	U4-5347-A	Rubber Foot	1
25	U4-1512-A	Terminal(+,-)	1
26	U4-5348	Exhaust Valve	1
27	08:4007416	Rating Label	1
28	CD-0.01UTZ	Capacitor	1
29	17:4006722	Washer	2
30	12:4006611	Manufacture Label	1
31	_____	Serial Label	1